


# CLICKIT ER PELVIS

## PELVIC FIXATOR

Universal clamp

Rapid application

Sterile kit

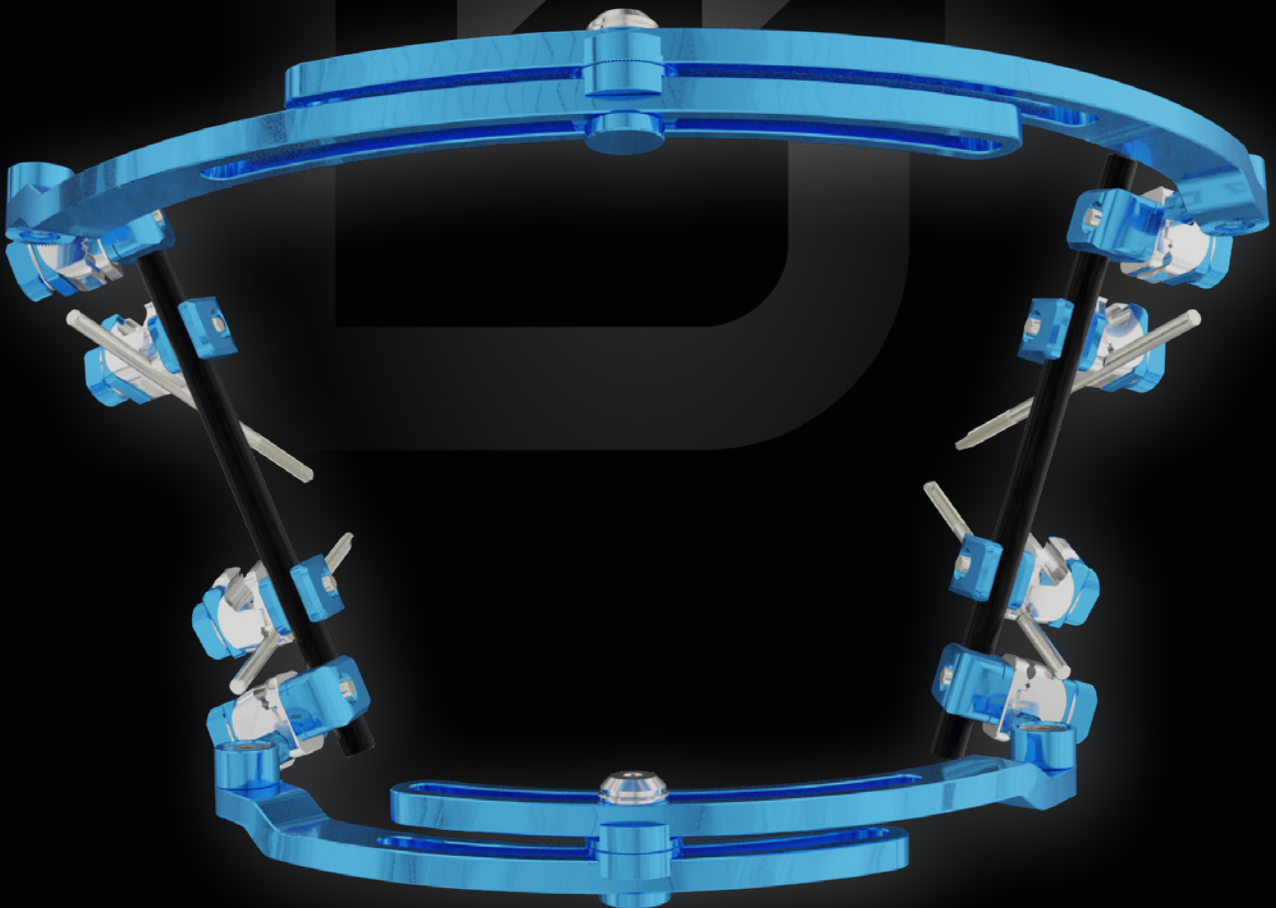

# Indications

Emergency stabilization of complete pelvic instability.

## FEATURES

## CLINICAL ADVANTAGES

Universal clamp

Multiplanar stabilization

Snap and click coupling

Sterile

rod/rod, rod screw (Ø 4.0-5.0-6.0 mm) and rod wire coupling

Increased stability

Quick application


Ready for use

# Components



### Universal clamp


One clamp for all types of coupling

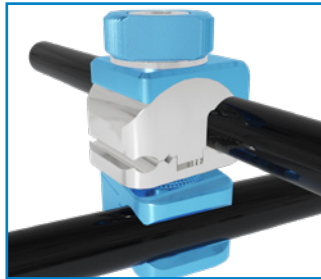

Rod/Rod Coupling



Rod/Screw coupling

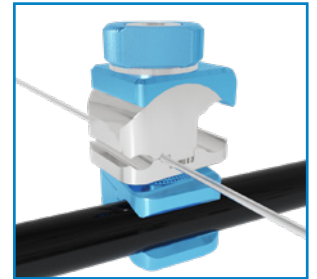

Rod/wire coupling




12mm Rods



### Pelvic arches

Made of aluminum alloy.  
Available in two sizes: large for the proximal side and small for the distal side

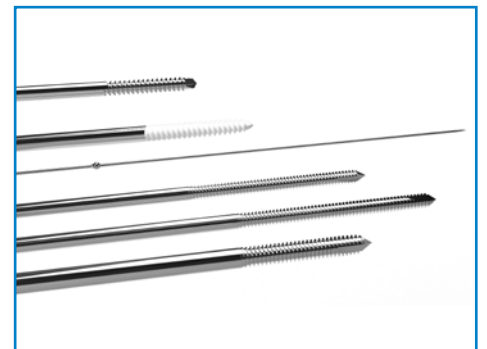

Self-drilling Cortical screws

# Kit ClickIt ER Pelvis



Description	Q.ty
Large arch (proximal)	1
Small arch (distal)	1
Universal clamps	8
Carbon rods Ø 12 mm x 250 mm	2
Self-drilling cortical screws Ø 6.0 mm L 250 mm thread 55 mm	4
Dedicated instrumentation	

## Configuration



## Sterilization

The components of the ClickIt ER system are supplied sterile.

We indicate sterilizing the instruments in saturated steam autoclave according to the following parameters:

**TEMPERATURE:** 134°

**EXPOSURE:** 5 min

# CLICKIT ER PELVIS 5006542ST

## COMPONENTS AND ACCESSORIES

Code	Description
5006508ST	Universal clamp
5006504ST	ER multiscrew clamp
5006510ST	Rod Ø 12 L 100 mm
5006512ST	Rod Ø 12 L 150 mm
5006514ST	Rod Ø 12 L 200 mm
5006516ST	Rod Ø 12 L 250 mm
5006518ST	Rod Ø 12 L 300 mm
5006520ST	Rod Ø 12 L 350 mm
5006522ST	Rod Ø 12 L 400 mm
5006526ST	Rod Ø 12 L 500 mm
5006528ST	Rod Ø 12 L 600 mm

## INSTRUMENTS \*

Code	Description	Q.ty
5006562	Complete ClickIt ER instrument trax	
Composition:		
5000124	3.0 mm hex screwdriver	1
5000222	Twist drill Ø 3.5 mm L 200 mm	1
5000600	T-wrench for screws Ø 4.0 mm	1
5000601	T-wrench for screws Ø 5.0 mm	1
5000602	T-wrench for screws Ø 6.0 mm	1
5000604	T-handle with AO attachment	1
5000605	Hexagonal rod Ø 5.0 mm AO connection	1
5006.04	Ø 4.0 mm adapter with AO connection	1
5006.05	Ø 5.0 mm adapter with AO connection	1
5006.06	Ø 6.0 mm adapter with AO connection	1
5000622	5.0 mm L-shaped Allen key	1
5006000	3.0 mm L-shaped Allen key	1
5006010	Drill guide Ø 3.5 mm L 110 mm	1
5006564	ER rod template	1
5006560ER	ClickIt ER sterilization trax	1

\* Implantable and sterile devices are covered by the CE-1936 marking, while non-sterile reusable instruments are covered by the CE marking





**HEADQUARTERS**

Via P. Gobetti, 56r - 16145 Genova  
Tel. +39 010 30801 - Fax +39 010 3080210  
[www.mikai.it](http://www.mikai.it) - [servizio.clienti@mikai.it](mailto:servizio.clienti@mikai.it)

**MANUFACTURING FACILITY**

Via Canestrello, 2 - 36050 Monteviale (VI)  
Tel. +39 0444 950100 - Fax +39 0444 950133  
[www.mikai.it](http://www.mikai.it) - [servizio.clienti@mikai.it](mailto:servizio.clienti@mikai.it)